

VitaPCR™


Real-time PCR results in 20 minutes



VitaPCR™

VitaPCR™ is an accurate and rapid molecular diagnostic solution that detects infectious diseases in just 20 minutes using real-time Polymerase Chain Reaction (real-time PCR) technology. It complements other laboratory-based methods and conventional lab testing, empowering healthcare professionals with highly accurate information at the Point-of-Care Testing (POCT).

Diagnostic Procedure in 4 Steps

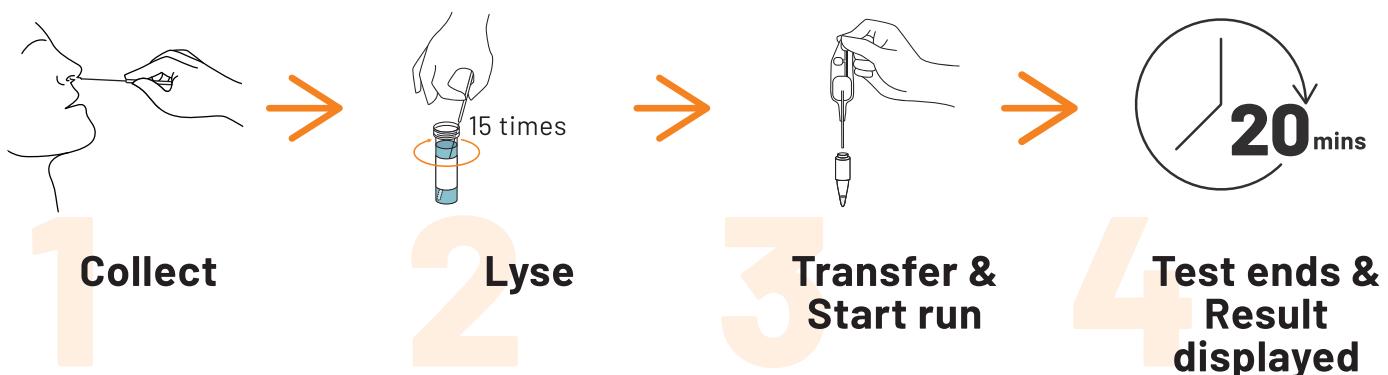

EASY TO USE

- Minimal training required
- Intuitive result interpretation
- Results with Ct value & fluorescence curve



RAPID

- Sample to result within 20 mins
- 1 min hands-on time for preparation
- Import patient ID via built-in scanner



VERSATILE

- Can be deployed anywhere
- No additional equipment needed
- Transportation and storage of reagents at room temperature



Hear from expert

The VitaPCR™ is a highly-sensitive test that delivers results in less than half an hour and to decentralize testing for SARS-CoV-2 in acute-care hospitals or emergency departments where rapid triage decisions are required.

Contribution of VitaPCR™ SARS-CoV-2 to the emergency diagnosis of COVID-19

Decentralized PCR Testing by Professionals



Emergency
Rooms



Nursing
Homes



Laboratories



Airports



Factories




Limited
Healthcare
Resources

VitaDataLink™

VitaDataLink™ is a data management software designed to connect to the VitaPCR™ platform. The software automatically transfers HL7 messages across a myriad of healthcare systems.

Intuitive Result Interpretation

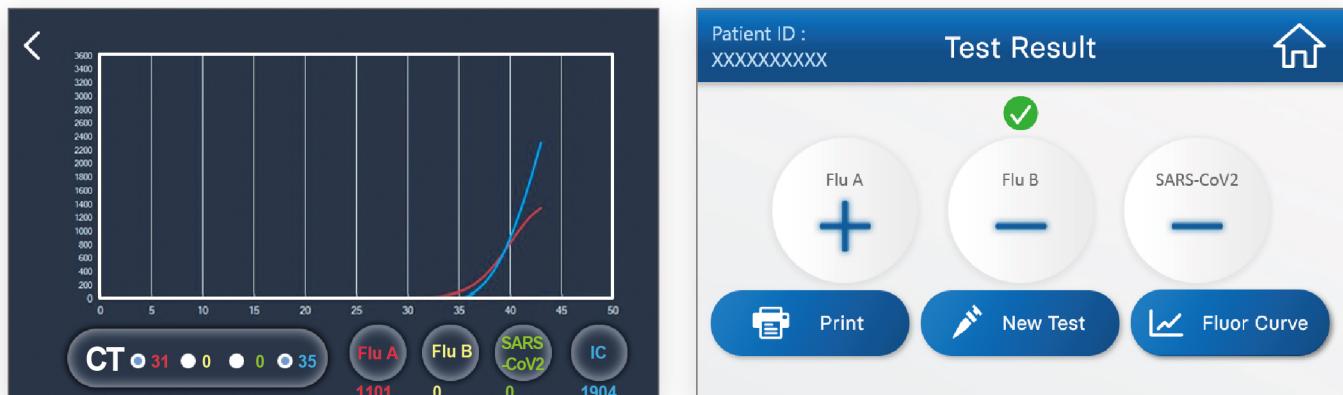

Diagnostic Menus

NS Nasal NP Nasopharyngeal OP Oropharyngeal

Assay	SARS-CoV-2 Gen2	Flu A&B	Flu A&B/ SARS-CoV-2	Flu A&B/RSV	Strep A
Catalog No.	PCRAE0120	PCRAE0103	PCRAE0128	PCRAE0106	PCRAE0107
Targets	N1, N2 genes	Influenza A: M gene Influenza B: NS gene	Influenza A: M gene Influenza B: NS gene SARS-CoV-2: N gene	Influenza A: M gene Influenza B: NS gene RSV: N gene	spy1258 gene
Sample Types	NP or OP swabs	NS or NP swabs	NP swabs	NS or NP swabs	OP or Throat swabs
Technology			Real-time PCR		
Hands-on Time			<1 minute		
Turnaround Time			20 minutes		
Kit Storage			5 - 25 °C / 41 - 77 °F		
Quantity			20 tests/box		

Detection of Variants

Credo Diagnostics Biomedical Pte. Ltd. continuously monitors the sequences of respiratory viruses listed in the major databases GISAID*, NCBI*, and BV-BRC* to provide confidence by tracking new emerging variants and broad subtypes. Based on BLASTn analysis comparing primers and probes of the assays respiratory viruses detection, the performances of the assays are not impacted by the variants and subtypes below.

SARS-CoV-2

VOIs

- XBB.1.5, XBB.1.16, FE.1.1, XBB.1.42.2
- EG.5, EG.5.1 (includes EG.5.1, HK.3 and HV.1)
- BA.2.86, JN.1(JN.1.11.1, KP1.1, KP.2, KP.3, LB.1, JN.1.16, JN.1.16.1, JN.1.7, JN.1.18)

VUMs

- XBB
- XBB.1.9.1, XBB.2.3

Previously Circulating

- Alpha variant: B.1.1.7
- Beta variant: B.1.351
- Delta variant: B.1.617.2
- Epsilon variant: B.1.427, B.1.429
- Eta variant: B.1.525
- Gamma variant: P.1
- Iota variant: B.1.526
- Kappa variant: B.1.617.1
- Lambda variant: C.37
- Mu variant: B.1.621
- Omicron subvariants: B.1.1.529, BA.1, BA.2M BA.3, BA.4, BA.5 (BQ.1.1 and BF.7), BA.2.75, CH.1.1, XBB.1.9.2, and DV.7
- Zeta variant: P.2

Influenza A

Subtypes


- H1N1, H1N2, H3N2, H5N1, H5N2, H3N8, H5N6, H7N3, H7N9, H9N2, H10N5

Variants

- H1N1 pdm09 variants¹
- H3N2 variants: seven variants belong to clade 3C.2a1b.2a.2 and one variant belongs to clade 3C.2a1b.1a as listed in the Danish published study²

50

Countries
Worldwide

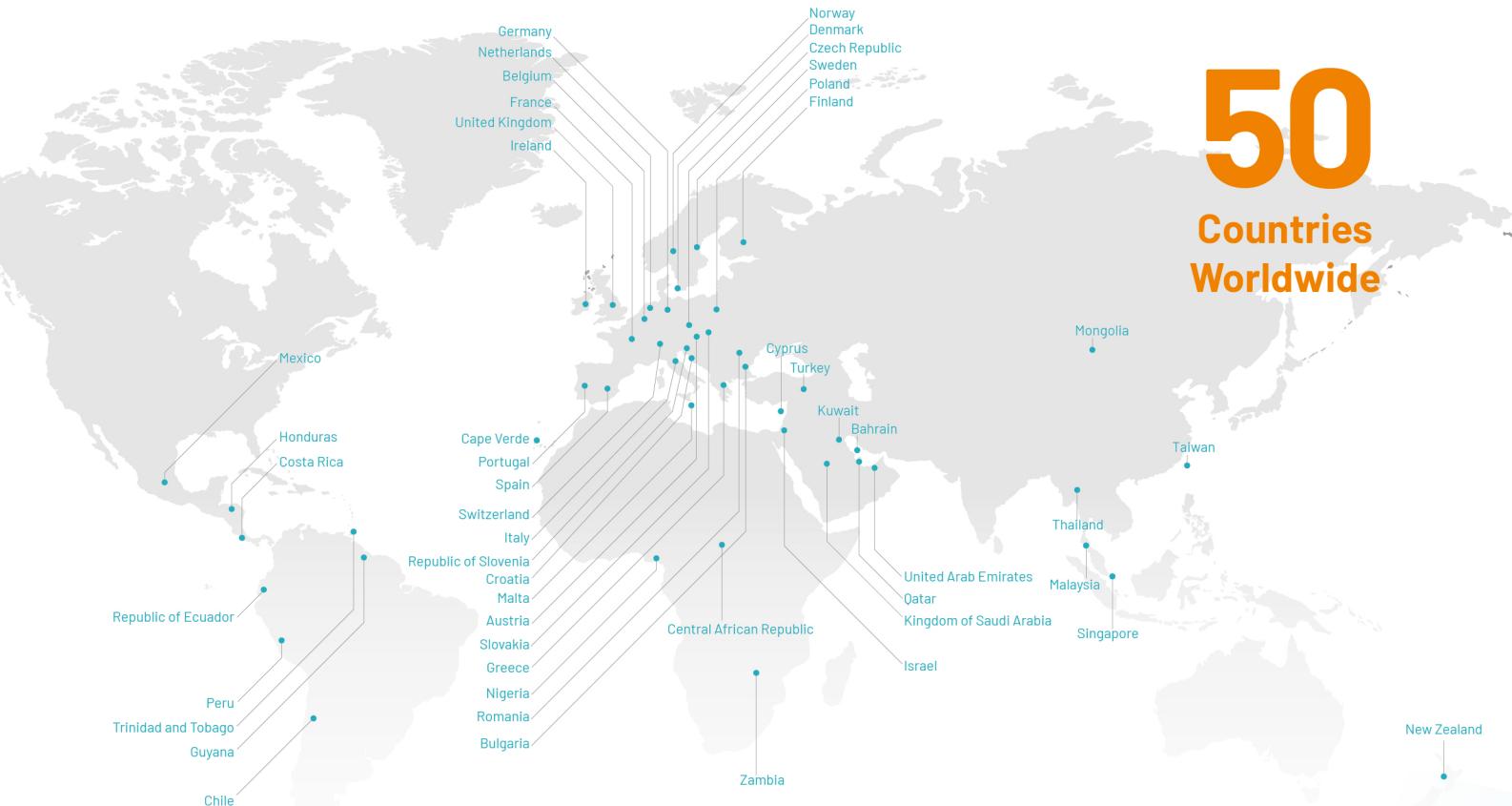

VitaPCR™ Instrument

Product Code	PCRAC0101
Dimensions	15.5 x 16.5 x 20.5 cm (H x W x D)
Weight	1.2 kg
Power Supply	12V, 5A
Color Touch Screen	4" capacitive touch LCD display
Fluorescence Detection	450 to 750 nm
Fluorescence Channel	3+1(Internal Control)
Storage Environment	15-40°C
Operating Environment	10-38°C, 10-80% RH
Power Adaptor	INPUT: AC 100-240V~2.0A Max, 50-60Hz OUTPUT: DC 12V-5A
Registration Status	CE-marked (IVDR)

Knowledge Empowering People

Credo Diagnostics Biomedical Pte. Ltd.
 491B River Valley Road #06-01 Valley Point Singapore 248373
[✉ service@credodxbiomed.com](mailto:service@credodxbiomed.com)
[🌐 credodxbiomed.com](http://credodxbiomed.com)
[in linkedin.com/company/credo-diagnostics](https://linkedin.com/company/credo-diagnostics)

cred
DIAGNOSTICS
Knowledge Empowering People